Cure Aqua®

Step 1: Prepare yourself and your environment

- Wash your hands with warm, soapy water, then dry them thoroughly.
- Get into a comfortable position near or on the toilet.
- Clean your genital area with warm, soapy water or a wipe.
- If necessary, gently pull back your foreskin to clean your penis, particularly around the entrance to your urethra.

Step 2: Prepare your product

- Use the dual finger loops on the top of the pack to peel it open easily.
- Use the sticky dot located on the back of the package to hang it vertically on a wall or smooth surface.
- Locate the protective no-touch sleeve. To keep the process hygienic and reduce contamination risk, only touch the sleeve and/or the colored funnel, do not touch the catheter itself.

Step 3: Insertion and voiding

- Hold your penis near the end and point it towards your stomach.
- Insert the catheter tip into your urethral opening.
- Slowly and steadily insert the catheter through your urethra.
- If you encounter slight resistance do not pull the catheter in and out, apply gentle but firm, continuous pressure until the catheter advances.
- While pointing your penis down, continue to push the catheter until urine begins to drain.
- Let the urine flow out of the catheter and into the toilet.

Step 4: Removal

- Once the urine has stopped flowing, you'll gently begin to remove the catheter from your urethra.
- If more urine begins to flow, stop moving the catheter until it finishes.
- You can then continue to remove the catheter, taking your time and ensuring you don't rush.
- You may notice a small amount of urine still comes out, but this is common and nothing to worry about.

Step 5: Disposal


- Once finished, remove the catheter and discard it immediately.
- Cure Aqua® is a single-use device and should be disposed of in the household waste.
- Be sure to use a new catheter each time you self-catheterize; do not reuse or reinsert the catheter.
- Thoroughly wash your hands when you're done.


convatec | CureAqua®


Male step-by-step user guide







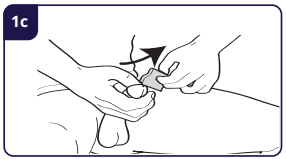

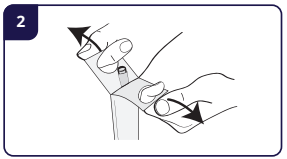


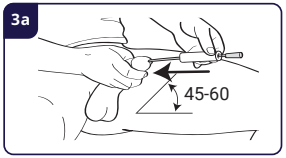


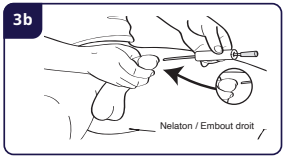

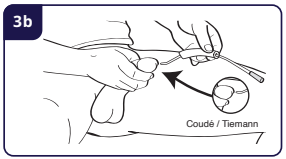



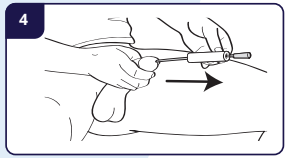

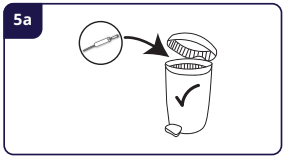



You're never alone with the Convatec me+ support program

Designed to help you live confidently with intermittent self-catheterization

convatec me+

What Convatec me+ offers you:

Educational tools & guides
A variety of resources to help you master your routine and feel confident.

Interactive video guides
Fully personalized experiences that show you how to catheterize and answer common questions.

Emotional well-being
A 10-module video program and workbook designed to help you navigate the emotional impact of catheterization.

Inspiration
Read stories from the Convatec me+ community of users who are living active, confident lives.

Speak with a member of the Convatec me+ support team today.
Call 1-800-422-8811 (Monday to Friday, 8:30 AM-7:00 PM ET).

Cure Aqua® is made of safe materials, as it does not contain potentially toxic or allergenic chemicals (no DEHP*/DINP, BPA, or natural rubber latex).[1]

*DEHP is used to soften PVC plastic and is suspected to reduce fertility. The chemicals are known to cause cancer, birth defects, and reproductive harm.[2]
BPA, bisphenol A; DEHP, di(2-ethylhexyl) phthalate; DINP, diisononyl phthalate; PVC, polyvinyl chloride.

REFERENCES

1. Data on file. Convatec. 2024. **2.** Proposition 65 Warnings. Di(2-ethylhexyl)phthalate (DEHP). https://www.p65warnings.ca.gov/fact-sheets/di2-ethylhexylphthalate-dehp. Accessed June 2026.

Convatec me+ is a trademark of Convatec.

AP-80926-GBL-ENGU-v2

convatec | CureAqua®

