
Wellspect®
A REAL DIFFERENCE

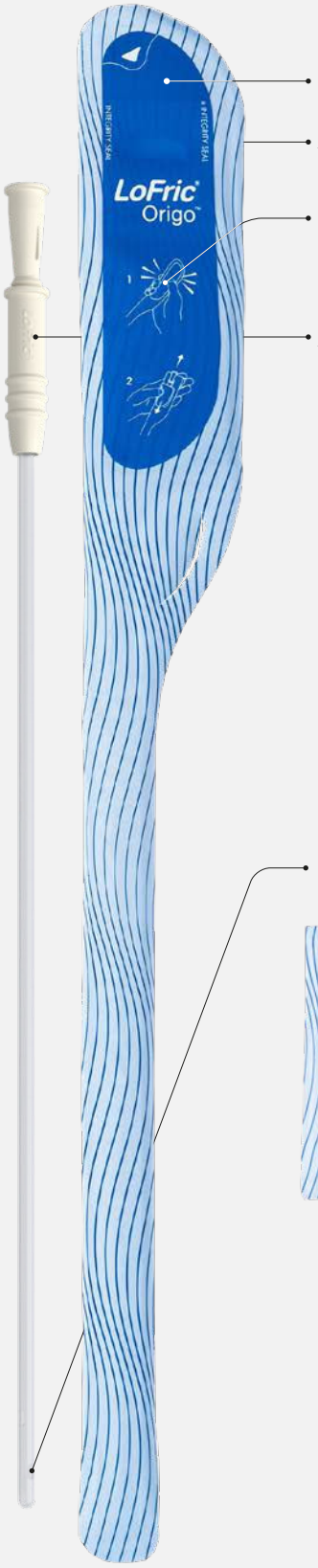

Peel-open tab.

Adhesive tab on back for reliable hanging.

Instant activation – sterile salt solution included, squeeze and it's ready to use.

Adjustable Insertion Grip allows better grip and control, promoting a non-touch technique.

Not manufactured with PVC, phthalates or latex.

Urotonic™ SURFACE TECHNOLOGY

LoFric catheters are coated with the unique Urotonic™ Surface Technology. This ensures the catheter remains smooth and safe, both during insertion and withdrawal.[1,2]

Smooth eyelets.[3,4]



Foldable into handy, discreet pocket size.



Nelaton *Coudé* *Flexible Coudé*

Product description

LoFric® Origo™, single use, hydrophilic urinary catheter with salt solution for activation.

Sterilization and shelf life

Sterilized by irradiation.
Use before expiry date on package.

Package

30 catheters per customer box. 4 customer boxes per shipping box.

Material composition

Hydrophilic surface coating: Polyvinyl pyrrolidone (PVP).

Catheter: Polyolefin-based elastomer (POBE).

Sachet: Laminate foil of polyethylene terephthalate (PET), aluminum and polyethylene (PE). Sterile water with sodium chloride (NaCl).

Item packaging:
Top foil: Laminate foil of polyethylene terephthalate (PET) and polyethylene (PE).

Bottom foil: Polyolefin.

Product range

Reference	Type	Length	Size
442 08 40	LoFric® Origo™ Nelaton (pediatric)	12”	08 FR
442 10 40	LoFric® Origo™ Nelaton (pediatric)	12”	10 FR
442 12 40	LoFric® Origo™ Nelaton (pediatric)	12”	12 FR
443 10 40	LoFric® Origo™ Nelaton (male)	16”	10 FR
443 12 40	LoFric® Origo™ Nelaton (male)	16”	12 FR
443 14 40	LoFric® Origo™ Nelaton (male)	16”	14 FR
443 16 40	LoFric® Origo™ Nelaton (male)	16”	16 FR
443 18 40	LoFric® Origo™ Nelaton (male)	16”	18 FR
444 10 40	LoFric® Origo™ Coudé (male)	16”	10 FR
444 12 40	LoFric® Origo™ Coudé (male)	16”	12 FR
444 14 40	LoFric® Origo™ Coudé (male)	16”	14 FR
444 16 40	LoFric® Origo™ Coudé (male)	16”	16 FR
444 18 40	LoFric® Origo™ Coudé (male)	16”	18 FR
447 12 40	LoFric® Origo™ Coudé FX (male)	16”	12 FR
447 14 40	LoFric® Origo™ Coudé FX (male)	16”	14 FR
447 16 40	LoFric® Origo™ Coudé FX (male)	16”	16 FR



User-friendly packaging.

NORDIC SWAN ECOLABEL 3098 0002

FSC


LoFric® Origo™

WELLSPECT does not waive any right to its trademarks by not using the symbols ® or ™.


Instructions

Single use hydrophilic urinary catheter with salt solution


Start by washing your hands with soap and water.

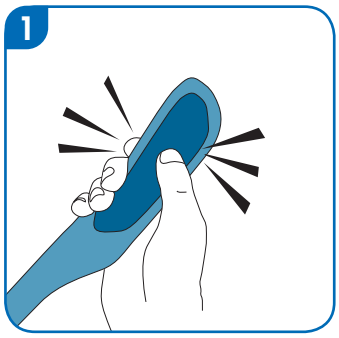


1. *Gently squeeze water pouch to release the salt solution down onto the catheter and the catheter is ready to use.*

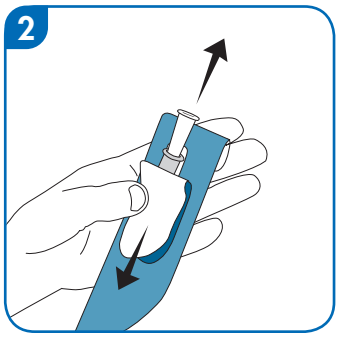

2. *Pull the tab down to open.*

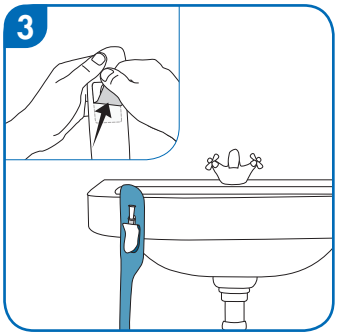


3. *Use the adhesive tab on the reverse side to hang the product, if necessary.*

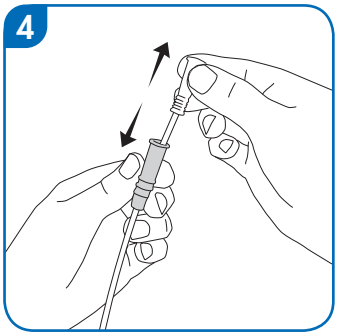

4. *Use Insertion Grip situated below the funnel for easy no-touch insertion.*
** Make sure to insert the catheter with the curved tip pointing upward.*

Discreet Carrying Option

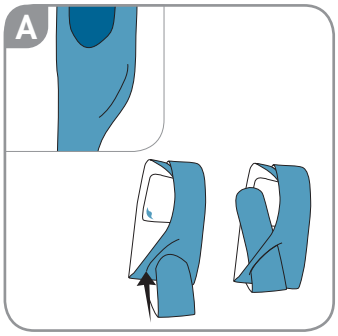

A. *Fold the package and use the slit to keep the catheter folded.*

IMPORTANT INFORMATION

This is a single use product. Once used, the surface coating will deteriorate and the catheter will no longer be sterile. Use only the salt solution included to wet the catheter. LoFric Origo must not be placed on paper or cloth as these materials can absorb the catheter's liquid coating.

© 2014 Wellspect HealthCare, a DENTSPLY International Company. All rights reserved. 75708-US-1301

Wellspect HealthCare, 840 Apollo Street, Suite 115, El Segundo, CA 90245
Phone: 1-877-4-LOFRIC (456-3742) www.lofric.us
marketinglofric-us@wellspect.com

Wellspect HEALTHCARE

Wellspect®
A REAL DIFFERENCE

References

1. Waller et al., Spinal Cord, 1997
2. Dunger et al., Acta Endocrinologica, 1988
3. Lundgren et al., Spinal Cord, 2000
4. Newman et al., Springer, 2018

At Wellspect we develop innovative continence care solutions that improve quality of life for people with bladder and bowel problems. We inspire our users to build self-confidence and independence as well as good health and well-being. We have been leading the industry for over 40 years with our product brands LoFric® and Navina™. We always aim to minimize the environmental impact of our products and passionately strive to become climate neutral. We work together with users and healthcare professionals to improve clinical outcome in a sustainable way, now and for the future.
Wellspect. A Real Difference.

For more information about our products and services, please visit Wellspect.us.
Join the conversation on LinkedIn, X, Facebook and Instagram.

wellspect.us

Wellspect HealthCare, 1235 Friendship Rd, Suite 205, Braselton, GA 30517.
Toll-free: +1 877 456 3742 .

CID68099-US-EN-2025-04-03 © 2025 WELLSPECT. All rights reserved.

LoFric

Manufacturer: Wellspect HealthCare, Aminogatan 1, P.O. Box 14, SE-431 21 Mölndal, Sweden.
Phone: +46 31 376 40 00. www.wellspect.com

CE 2797

℞ Only